

Adenoviral Knockdown Virus Data Sheet

Vector Information

Vector: E1/E3 deleted adenoviral serotype 5 vector
Expression cassette: shRNA under control of human U6 promoter

Target Information


Target Gene: **Caspase3**
Specimen: Sus scrofa
RefSeq#: NM_214131

Validation

Validation Method: Western Blot
Cell line used for validation: IPI21
Total protein prepared: 60 h post transduction
Antibody used for validation: CST antibody #9662
Detection Method: ECL
Quantification Method: Densitometrically
Quantification Software: Quantity One (BioRad)

Results

KD rel. to none transduced cells 70%

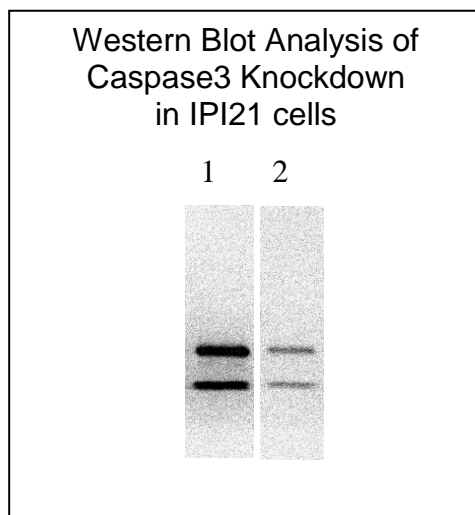

Western Blot showing Caspase 3 expression in extracts prepared from IPI21 cells (lane 1) and with the adenoviral vector Ad-shCaspase3-SS2 expressing a shRNA directed against Caspase3 mRNA (lane 2).