

Adenoviral Knockdown Virus Data Sheet

Vector Information

Vector: E1/E3 deleted adenoviral serotype 5 vector
Expression cassette: shRNA under control of human U6 promoter

Target Information


Target Gene: **Caspase6**
Specimen: Homo sapiens
RefSeq#: NM_001226

Validation

Validation Method: Western Blot
Cell line used for validation: A549
Total protein prepared: 60 h post transduction
Antibody used for validation: CST antibody #9762
Detection Method: ECL
Quantification Method: Densitometrically
Quantification Software: Quantity One (BioRad)

Results

KD rel. to none transduced cells 80%

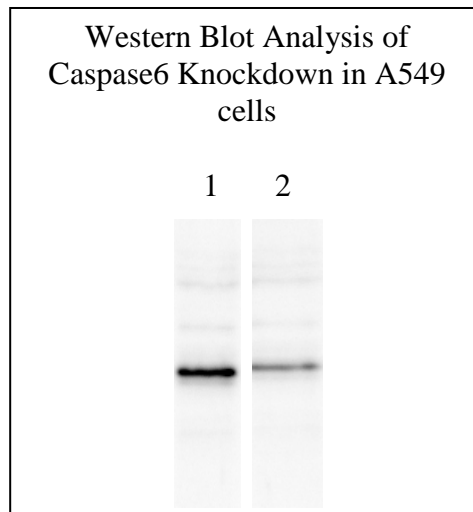

Western Blot showing Caspase6 expression in extracts prepared from A549 cells (lane 1) and with the adenoviral vector Ad-shCaspase6 expressing a shRNA directed against Caspase6 mRNA (lane 2).